

## HANDLING

1. Do not use liquid chemicals to clean the product, because it can damage the plastic parts of the lamp
2. To avoid short circuits and damage to the product, do not open the lamp in addition to technicians



## HALOGEN MEDICAL EXAMINATION LIGHT

**MODEL : MEL-01**



**Manual Book**

## CAUTIONS

1. Make sure that the resource has the appropriate voltage 220V-230V before it is connected to a power source
2. Keep away from water or damp place for proper use

## HOW TO USE

1. Plug the AC cord to the outlet with a working voltage of 220V

2. Press the ON button to turn on the light switch potensio



3. Turn pot light settings halogen lamp (arrangements of minimal, dim up to maximum)




4. Press the OFF button and remove the power cord from the wall socket when not in use halogen bulbs

## SPESEIFICATIONS

1. Voltage : 220 V
2. Power : 50 W
3. frequency : 50Hz / 60 Hz
4. Illuminance :  $\pm 44.000$  lx

## PARTS NAME

1. Button ON/OFF
2. Clinch Halogen Lamps
3. Cover Halogen Lamps
4. Circuit Casing
5. Circuit Cover Casing (to link and lock between the upper and lower lamp post)
6. Lower Lampposts
7. Wheel

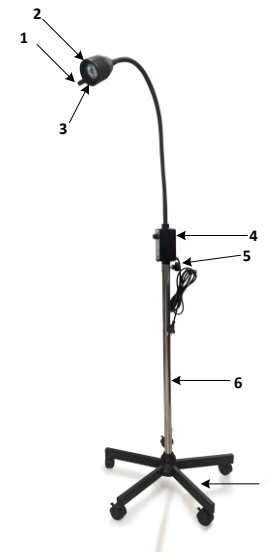

CONTROLLED COPY