

This user manual will help you understand the operation and maintenance of the product properly. Please note that the product must be used correctly and must comply with this manual. Users who do not comply with this manual may result in the unit not functioning properly or an accident in which PT. Sinko Prima Alloy cannot be held liable. PT. Sinko Prima Alloy owns the copyright of this manual. Without prior written approval from PT. Sinko Prima Alloy. All materials related to this manual may not be photocopied, reproduced, or translated into other languages.

Materials are protected by copyright law, but not limited to confidential information such as technical information and patent information contained in this manual, users may not provide information to third parties that are irrelevant.

The user should understand that nothing in this manual provides express or implicit permission to the user. All rights or licenses to the intellectual property of PT. Sinko Prima Alloy.

PT. Sinko Prima Alloy reserves the right to modify, update and explain this manual.

#### Manufacturer's responsibility

PT. Sinko Prima Alloy only considers itself responsible for any effect on the safety, constraints and equipment performance of the unit, if:

The assembly, extension, re-adjustment, modification or repair operations are carried out by a person authorized by PT. Sinko Prima Alloy. The instrument is used according to the instructions for use. PT. Sinko Prima Alloy will provide a request for a part list, description, calibration instructions or other information on individual repairs to repair components of the unit appointed by PT. Sinko Prima Alloy as a component that can be repaired by service personnel.

#### Product Information

**Product name** : Elitech Baby Length Measuring Board /  
Papan Pengukur Panjang Bayi

**Model** : MTR BABY-002

#### Term in Condition

This manual is intended to provide a basic concept of security measures.

#### Warning

**A warning sign** suggest certain actions or situations that result in personal injury or death.

#### Attention

**An attention sign** suggest actions or situations that result in product damage, inaccurate data or cancellation of procedures.

#### Note

**A note** prove useful information in terms of functions or procedures.

## Introduction

### The intended use


MTR Baby-002 is a device that is intended to determine the length or height of a baby, the range that can be measured is 10 ~ 100 (cm).

This tool is designed to be very simple, because in its use, the baby can simply be put to sleep on the device with the head attached to the headbase and then place the footpiece into the baby's feet. This tool is suitable for babies in the age range of 0 ~ 23 months.

### Feature:

- Easy operation.
- Resolution 1 mm
- Accurate reading

Part :

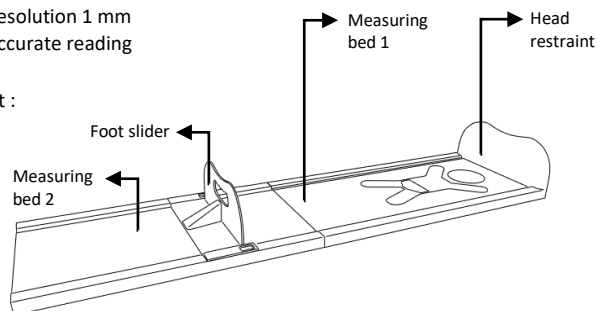

Pict 1

## Preparation and How to Use

1. Arrange measuring bed 1 and 2, insert them into the locking knob
2. Position the baby in a supine position with the head on the head restraint. Position baby's gaze at 90° up
3. Hold the baby's knees / do not let the baby's legs bend.
4. Slide the foot slider until it touches the baby's feet.
5. Measurement can see the scale shown.

## Maintenance and Care

1. After use, please clean the baby's measuring board.
2. Store the baby's measuring board in a closed place, at room temperature so that the device does not expand.
3. Avoid stacking measuring tools.
4. Avoid hard impact.
5. Periodic maintenance.

## Attention

1. Provide service to a professional technician.
2. Store the device in a clean environment and avoid vibration during storage.
3. Do not clean the surface of the unit with liquids that can cause corrosion.
4. Every use must be cleaned, can use medical alcohol and dry cloth.
5. Place the unit onto a flat media and minimize the use of the base.
6. If using another tool or base, adjust the height of the unit again so that it is equal to or equal to the zero point.
7. The position or posture of the patient's body is very influential on the reading, read the instructor's use again for measurements.

## Specification

Model	MTR BABY-002	
Range	10 ~ 100 (cm)	
Material	Body	ABS
	Foot pad	Rubber
Resolution	1 mm	
Dimension	1120 mm(L) x 300 mm(W) x 150 mm (H)	
Weight	≤ 2,5 kg	
Operating condition	Temperature	10°C ~ 40°C
	Humidity	30% ~ 80% RH
	Atmosphere pressure	800 hpa ~ 1050 hpa
Transport and storage	Temperature	-20°C ~ 40°C
	Humidity	≤95% RH
	Atmosphere pressure	800 hpa ~ 1050 hpa

## Warranty and Service

### Warranty

PT. SINKO PRIMA ALLOY guarantees that PT. SINKO PRIMA ALLOY has complied with the table specifications of the product and will be free from defects in materials and workmanship that occur within the warranty period.

Warranty does not apply if:

- Damage occurs due to handling when shipping expeditions.
- Damage caused by improper use or maintenance.
- Damage caused by changes or repairs by anyone not authorized by PT. SINKO PRIMA ALLOY.
- Damage caused by accident.
- Replacement or removal of serial number labels and manufacture labels.


If a product covered by this warranty is determined to be defective due to faulty materials, defective components or workmanship and a warranty claim is made within the warranty period, PT. SINKO PRIMA ALLOY at its sole discretion repairs or replaces damaged parts free of charge. PT. SINKO PRIMA ALLOY will not provide a replacement product for use while the defective product is being repaired.

### Contact Information


If you have any questions regarding maintenance, technical specifications or device malfunction, please contact your local distributor.

## Reference


Length-for-age GIRLS  
Birth to 6 months (z-scores)

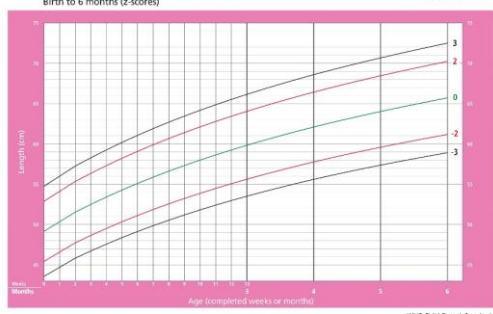


Length-for-age GIRLS  
6 months to 2 years (z-scores)

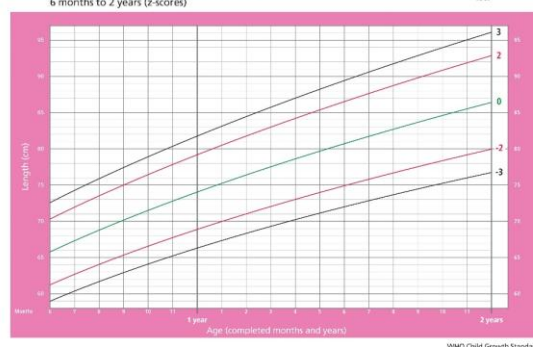

Length-for-age BOYS  
Birth to 6 months (z-scores)

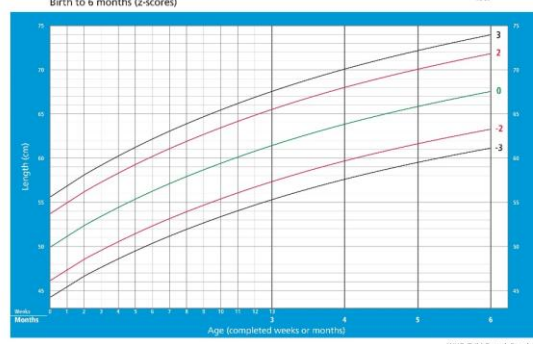

Length-for-age BOYS  
6 months to 2 years (z-scores)

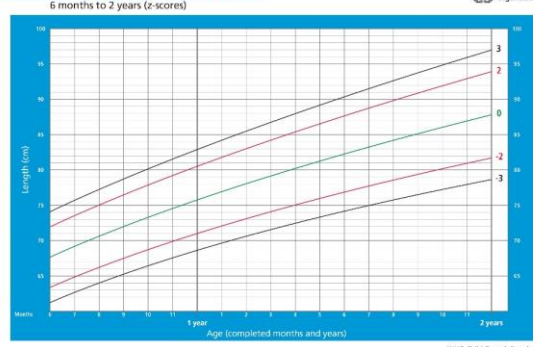

### PT. SINKO PRIMA ALLOY

TAMBAK OSOWILANGUN NO.61  
PERGUDANGAN OSOWILANGUN PERMAI BLOK E7-E8 SURABAYA - 60191

No. Telephone 031-7492882,74828816,7482835

Email: sinkoprima@gmail.com

teknik.sinkoprima@gmail.com

Website: www.elitech.id

Nomor Dokumen : SPA-BM/PROD-165

Tanggal terbit : 23 Mei 2023

Rev : 05