

PT. SINKO PRIMA ALLOY

TAMBAK OSOWILANGUN NO.61
PERGUDANGAN OSOWILANGUN PERMAI BLOK E7-E8
SURABAYA - 60191

TLP. 031-7492882,74828816,7482835

sinkoprima@gmail.com

teknik.sinkoprima@gmail.com

Website: <http://www.elitech.co.id>

Nomor Dokumen : SPA-BM/PROD-66

Tanggal Terbit : 19.11.2019

Revisi : 01

File Ver : 1.02

TIMBANGAN BAYI MEKANIK

BABY DIGIT-ONE

CONTROLLED COPY

MANUAL BOOK

Thank you very much for buying a product from us Elitech Mechanical Baby Scales (Baby Digit-one). Please read the following manual before operating this instrument.

For Users

© Copyright by PT. Sinko Prima Alloy

PT. Sinko Prima Alloy

Industrial complex, Jl. Osowilangun Pond No. 61

Osowilangun Permai Warehouse Block E-8.

Surabaya -East Java. Indonesia

Tel: + 62-31-7482816, 7492882 Fax: + 62-31-7482865

Postal code: 60191

E-mail: sinkoprima@gmail.com

Website: <http://www.elitech.co.id>

CONTROLLED COPY

5. Attention safety

- The scale should be placed on a hard, flat surface for the baby's comfort and safety.
- Keep the scale clean to ensure normal work.
- Place the baby in the center of the tray, otherwise it might affect the weight accuracy.
- Overload (heavier than 25kg) is strictly prohibited.

Note

- ✓ PT. Sinko Prima Alloy product properly, please read the User Manual before using this product.
- ✓ Please save the User Manual carefully.
- ✓ The company is not responsible or does not provide free maintenance for any faults or damages because the user does not follow the instructions from the manual and does not maintain and maintain the equipment.
- ✓ This company provides quality products according to company standards.

A Warning

- avoid exposure to hard objects, moisture, dust, chemicals, liquids, large temperature fluctuations and proximity to heat sources (open flames, radiators).
- Repairs can only be done by PT. SINKO PRIMA ALLOY
- Cleaning: You can clean the scale with a damp cloth and a little liquid, if needed. Do not immerse the scale in water or rinse it under running water.
- Storage: Do not place anything on the scale when not in use.
- Strong electromagnetic fields can interfere with the accuracy of your balance.
- This scale can accept a maximum load of 25Kg.
- Placing the scale on a flat, solid surface is important to ensure that the scale can operate properly.

Operation Manual

1. Features

- Professional body scales, specially designed for health facilities
- Material per carbon manganese high quality, high stretch and accuracy.
- The 170mm diameter round and large needle display makes it easy to read.
- Measuring capacity up to 25kg.
- Graduation 50/100 grams.
- Commercial grade (accuracy class M3).
- Pure and food grade PP material tray.
- Sturdy and safe for babies.
- The curved tray fits the baby's anatomy and is comfortable.
- Non-slip and flexible rubber.


1. Spesifikasi teknis

- Model : Baby Digit - One
- Measuring Ability: 0gr - 25Kg
- Scale capacity: 25kg.
- Reading gradient: 50gr.
- Scale unit dimensions: 542mm (L) X 436mm (W) X 140mm (H)
- Tray height : 80 mm
- Weight : 3,35Kg

3. Parts of the scale

PT. Sinko Prima Alloy

02

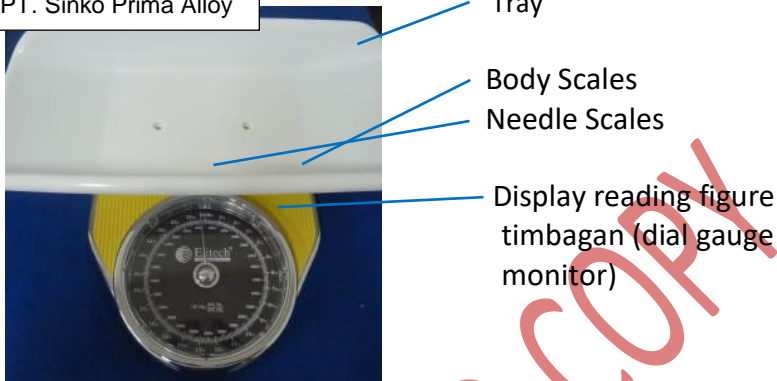

4. Intruksi Pengoperasian

- Install the tray / tray on the body of the scale, make sure the position of the screw tray hole is right on the screw hole in the body of the scale.
- Install two screws on the tray to lock the tray on the body of the scale, after tightly closing the screw with the rubber cover that has been provided.
- Make sure the scale needle shows zero before weighing. If the scale needle is not right at zero then turn the scroll on the back of the scale slowly to the right or left to position the scale needle to zero.

Caution: The scroll loop must not exceed one complete loop.

- Place the baby on the tray / tray and the reading of the scale is shown on the display.

CONTROLLED COPY